Passport  
*for the «ImperStar» device*

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## Contents

### EN

1. General guidance.....	4
2. The purpose of the device.....	4
3. Technical data.....	4
4. Completeness of delivery.....	4
5. Device.....	5
6. Preparation for work.....	5
7. How to use the device.....	6
8. Precautions when working with the device.....	7
9. Maintenance of the unit.....	8
10. Typical malfunctions and methods of their elimination.....	8
11. The device acceptance card.....	9
12. Documents of the Ministry of Health of Ukraine.....	10
13. The passport for the « <i>ImperStar</i> » device.....	11

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## 1. General guidance

This operating instruction covers the «*ImperStar*» device (hereinafter - the device) which is intended to familiarize with the purpose of the device, its design, operation, study the rules of operation and maintenance of the device.

## 2. The purpose of the device

2.1. The device is designed for:

- active bioresonance correction based on the generation of weak electromagnetic waves of a special shape and frequency, which are transmitted through the body with the help of the device;
- 2.2. The effect is carried out by waves having a biomagnetic nature with a corrective soft, non-depressing effect on the nucleus and plasma of the cell and its components.
- 2.3. The action of the device aimed at correcting the homeostasis of the human body is one of the few effective ways to prevent the disease.
- 2.4. The device is intended for use by specialists of beauty salons, massage therapists, rehabilitation specialists, etc., as well as for self-use by consumers at home.

## 3. Technical data

- 3.1. According to DSTU 3798 the device refers to the products with internal power supply (type of protection against electric shock is B).
- 3.2. The device operates autonomously without connection to the power supply network.
- Built in lithium polymer voltage battery 3.7V is used as power supply.
- 3.3. During operation the device should be resistant to climatic factors in accordance with GOST 20790 for climatic performance UHL 4.2:
- at ambient temperatures from 10°C to 35°C;
  - at a relative humidity of 80% at a temperature of 25°C.
- 3.4. Current consumption not exceeding 30 mA at 3,7 V supply voltage.
- 3.5. Mean time between failures is not less than 1000 hours.
- 3.6. The full period of service is not less than 5 years.
- 3.7. Operating mode setting time is not more than 5 seconds.
- 3.8. Continuous operation time is not less than 5 hours.
- 3.9. Overall dimensions, up to, mm — 105x65x25.
- 3.10. Frequency range 150 Hz to 400 kHz.
- 3.11. Weight without packaging - up to 0.4 kg.

## 4. Completeness of delivery


- 4.1. «*ImperStar*» device — 1 pc.
- 4.2. The power adapter — 1 pc.

- 4.3. USB cable – micro USB— 1 pc.
- 4.4. Passport - 1 pc.
- 4.5. Brochure on the use of the device - 1 pc
- 4.6. Packing container - 1 pc
- 4.7. Warranty - 1 pc
- 4.8. Set of flyers - 1 set.

## 5. Device

### 5.1. Controls and indications:

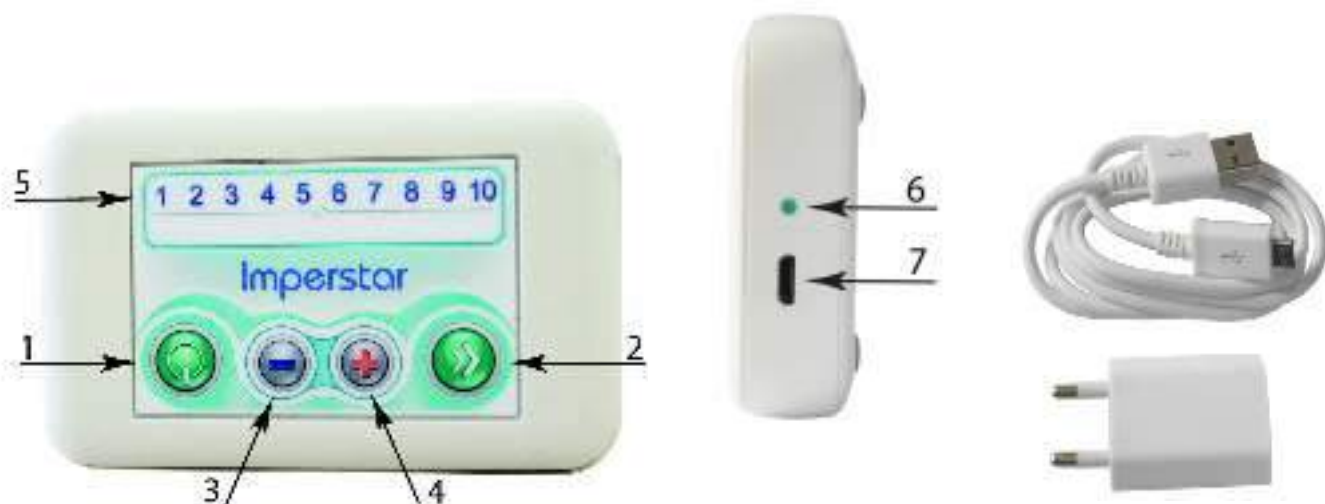- 1. System power on/off button.
- 2. Mode (program) selection button.
- 3. Button for lowering the signal strength.
- 4. Button to increase the signal strength.
- 5. Indicators of the selected mode (program).
- 6. Charge indicator) red - "charge", green - "charge end"
- 7. Micro USB connector

## 6. Preparation for work

- 6.1. After transporting or storing the device at low temperatures, the device should be kept at room temperature for at least 2 hours before switching on.
- 6.2. Remove the device from the packaging and make sure by external inspection that there is no mechanical damage. Check the completeness of the device.

## 7. How to use the device

- 7.1. Press the power on/off button (1) three times to turn on the device. The program selection indicator (5) should light up (in running mode). During the operation, the device beeps every minute of the session.
- 7.2. Use the program selection button (2) to select the desired program. In this case, the corresponding program indicator will blink and light up with a fixed light, which will send out a double sound signal

### 7.3. Process 1, 2, 3, 4, 5.

- 7.3.1. Turn on the equipment according to 7.1.
- 7.3.2. According to paragraph 7.2, select the appropriate procedure from 1 to 5.
- 7.3.3. Put it in your body or in your pocket.
- 7.3.4. The working time of the programme is seven minutes. The device sends out sound signals every minute. After 7 minutes, the device enters a 20 minute rest state.
- 7.3.5. After 20 minutes, the device will start to send out the sound signal, and then can continue the conversation by pressing the program selection button once) 2. The whole conference cycle: 7 minutes work, 20 minutes rest, 7 minutes work, 20 minutes rest, 7 minutes work, 7 minutes work.

### 7.4. Processing procedure № 6,7,8,9,10

- 7.4.1. Turn on the equipment according to 7.1.
- 7.4.2. According to paragraph 7.2, select the appropriate procedure from 6 to 10.
- 7.4.3. Put it in your body or in your pocket.
- 7.4.4. Massage with soft motion of clockwise rotation (
- 7.4.5. Working hours range from 7 minutes to 28 minutes), 1 to 5 times a day if necessary.

### 7.5. Transfer of information to water.

In order to increase the effective and mild effects on the body, it is recommended to transfer information from the device to the water.

This device can be used to record information on both a small volume of liquid and a larger volume.

#### Using a small container (maximum liquid volume is 1 liter)

This requires:

- 7.5.1. Place the device on a flat surface with the front panel facing up.
- 7.5.2. Turn on the device and select the program you want to record on the water (the signal strength of "+" or "-" signal is not necessary to be chosen)
- 7.5.3. Install a small container of liquid (maximum volume is up to 1 liter) on the front panel.

**Attention!** Recording the program on the water takes 7 minutes, after which the liquid is ready for use.

**Attention!** A common consumer mistake. If you want to record several programs in a row, only the last program is recorded.

**Use of large volume (1 liter or more)**

This requires:

7.5.4. Place the device on a flat surface, set to the desired volume, with the front panel facing up.

7.5.5. Turn on the device and select the program you want to record on the water (the signal strength of "+" or " - " signal is not necessary to be chosen)

*Attention! To avoid failure of the device, do not immerse the device in liquid.*

**TEN PROGRAMS OF THE DEVICE**

№	Name of programs
1	Antifungal.
2	Antibacterial and antiviral.
3	Anthelmintic and Antiprotozoal (against protozoa).
4	General, including anthelmintic, antiviral (including Human Papilloma Virus and verruga), antiprotozoal (against protozoa), and miticidal.
5	Preventive and detoxification.
6	Cosmetology and trophism improvement.
7	Analgesic and anti-inflammatory.
8	Rehabilitation and prevention of the musculoskeletal system.
9	Improvement of blood vessels, anti-cellulite, lymphatic drainage.
10	Improvement and prevention of problems of the gastrointestinal tract.

**8. Precautions when working with the device**

8.1. According to DSTU 3798 the device refers to the products with internal power supply (type of protection against electric shock is IIa).

8.2. The device is electrically safe, no grounding is required.

**It is strictly forbidden:**

- to perform any procedures using a faulty device;
- to carry out repairs when the device is switched on.

8.3. Working with the device is allowed only after reading the operating Instructions and electrical safety rules.

8.4. In case of malfunction — switch off the device immediately by pressing the on/off button (1).

## 9. Maintenance of the unit

9.1. Maintenance and testing of the device operability should be performed by the consumer.

9.2. It is required to protect the device from moisture and aggressive environments, protect from shocks, and sudden temperature changes.

9.3. Keep the unit away from direct sunlight.

**Attention!** To prevent damage to the device, do not use solvents or abrasives to clean the outside of the device.

## 10. Typical malfunctions and Troubleshooting

№	Name of a malfunction and its possible cause external manifestation	Possible reason	Troubleshooting method
1.	When switching on the device does not light up There are no indication elements. The power supply unit does not switch on	Battery discharge	Charge the battery
2.	The button did not respond, there was no specific click	Button failure	Contact Service Center
3	When you turn on the instrument and program, turn on "run light mode," and none of the lights are on fire.	Led failure(Ability to work is not affected)	Contact Service Center

## 11. The device acceptance card

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№\_119\_\_\_\_\_

Manufactured and accepted according to the mandatory standard requirements, acting normative instruments and is classified as fit for use.

The system was manufactured on the basis of TU U 26.2-43746219-001:2020  
Manufacturer: «Impermaster» LTD, Ukraine, Kharkov, Moskovsky ave. 273

Chief of QCD \_\_\_\_\_ ( \_\_\_\_\_ )  
(personal signature)

P. P.

*Dated of acceptance* « \_\_\_\_\_ » \_\_\_\_\_ 20 \_\_\_\_\_

## 12. Documents of the Ministry of Health of Ukraine

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